

## ARCADIS Orbic

**SP**

### Installation Instructions

System

Unpacking the Basic Unit

© Siemens AG 2004

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

**Document revision level**

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

**Disclaimer**

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

1

Unpacking the Basic Unit

4

Work aids . . . . . 4

workflow . . . . . 5

Changes to previous version . . . . . 9

## Work aids

Required tools:

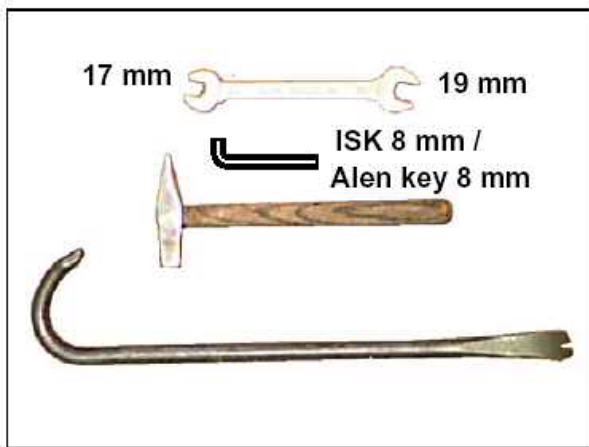


Fig. 1: Tools\_

## workflow

- Cut the fastening strap and remove the ramp from the pallet ([Fig. 2 / p. 5](#)):



Fig. 2: Attachment

- Unscrew and remove the retaining clamp ([1/ Fig. 3 / p. 5](#)).

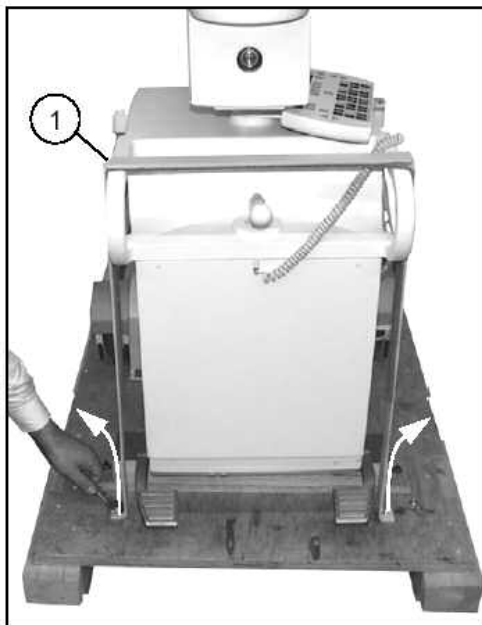


Fig. 3: Retaining clamp\_1

- Unscrew the holding block (Fig. 4 / p. 6).

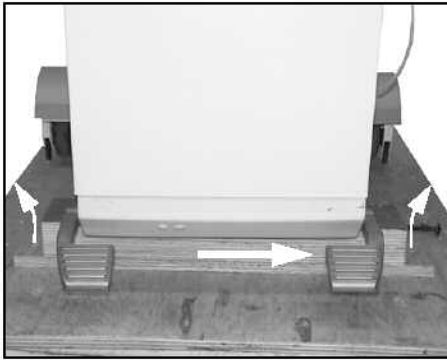


Fig. 4: Holding block

- Unscrew the retaining clamp (2/Fig. 5 / p. 6) on the left and right (two nuts on each side)

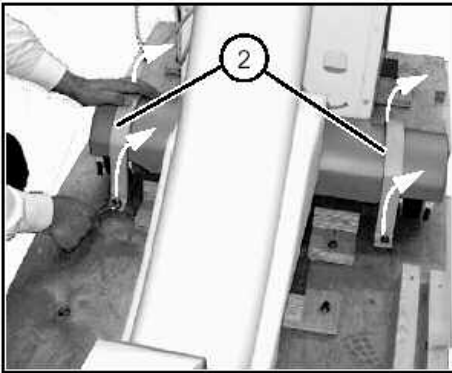


Fig. 5: Retaining clamp\_2

- To remove the transport clamp (1/Fig. 6 / p. 6), the front fastening screws (2/Fig. 6 / p. 6) of the railings must be **loosened** on both sides of the C-arm.
- Unscrew and remove the transport clamp (1/Fig. 6 / p. 6).
- Subsequently retighten the fastening screws (2/Fig. 6 / p. 6).

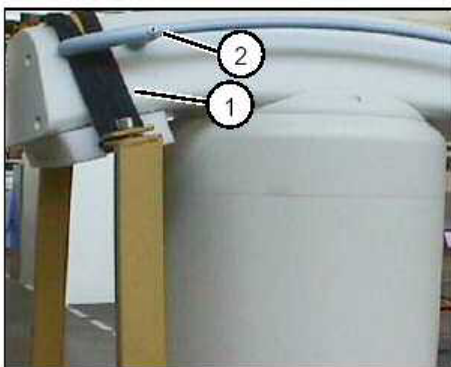


Fig. 6: Transport clamp

- Loosen and remove the support (Fig. 7 / p. 7).



Fig. 7: Support

- Unscrew and remove the retaining clamp (3/Fig. 8 / p. 7) and holding block at the tube.
- Unscrew, tilt, and remove the retaining clamp (4/Fig. 8 / p. 7).

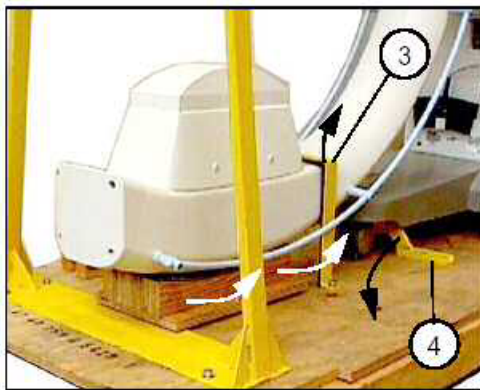


Fig. 8: Retaining clamp\_3\_4

- Position the ramp against the pallet and nail it down (Fig. 9 / p. 7).

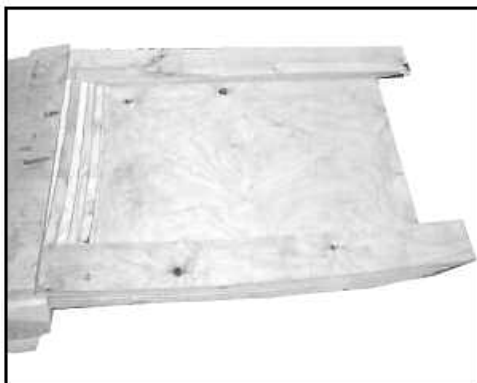


Fig. 9: Ramp\_pallet

- Release the foot brake and pull the unit in the direction of the pallet ramp (Fig. 10 / p. 8).



*Fig. 10: Foot brake*



## Changes to previous version

n.a.; initial version

